

Starlight® Protect-Evo



The Starlight® Protect-Evo builds on our experience with the Starlight® Protect and Protect Plus. The back-of-the-head area is extensively equipped with double-glued high-tech foam. The foam is ergonomically adapted to the shape of the back of the head. The preformed back-of-the-head area provides even more protection, especially for users who are at risk from uncontrolled fall directions. Impact forces are optimally absorbed and distributed.

The round shape of the ear straps guarantees an optimal fit of the chinstrap. The chinstrap can be easily adapted to the head shape of the user due to its free mobility on the ear strap. This guarantees the perfect fit of the head protection.

The protective effect in the upper head area is achieved by longitudinal and transverse braces made of shock-absorbing high-tech foam.

The thickness of the forehead protection is 1.5 cm.

Choose from various leather colours – our high quality standards apply also for our coloured leather!

The calfskin used for the internal lining is perforated. Thereby it is breathable, regulates heat and humidity and is very skin-friendly.



Head protection Starlight® Protect-Evo. The round shape of the ear straps guarantees an optimal fit of the chinstrap.

Technical details

Starlight® Protect-Evo	
Forehead protection	1.5 cm
Head circumference*	Made to measure
Padding	Hightech foam
Outer material	Leather
Lining	Perforated calf leather
Product weight (basic model)	approx. 230 g

*For technical reason, producing head protection in sizes under 48 cm in circumference is not possible.



Starlight® Protect-Evo – the back-of-the-head area is extensively equipped with double-glued high-tech foam adapted to the shape of the head.



*Starlight® Protect-Evo with ear protection and
alternativ upper side*




Starlight® Protect-Evo with ear protectors – side view

Models / Accessories

Starlight® Protect-Evo
Basic model

Article number: 5782-...

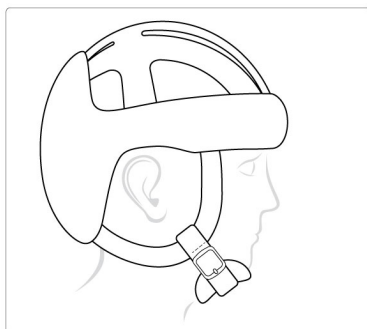


Starlight® Protect-Evo
Basic model

Article number: 5782-...

...with alternative top

Article number: 5789-...

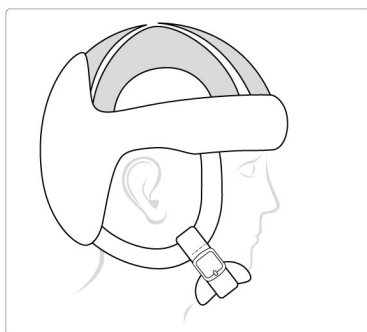


...with alternative top

Article number: 5789-...

...with ear protection

Article number: 5784-...

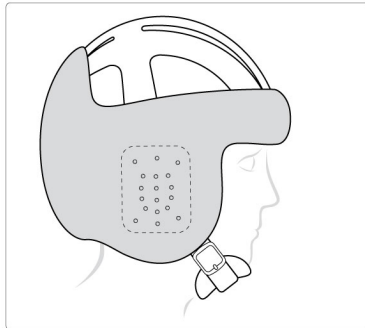


...with ear protection

Article number: 5784-...

**...with alternative top
and ear protection**

Article number: 5791-...

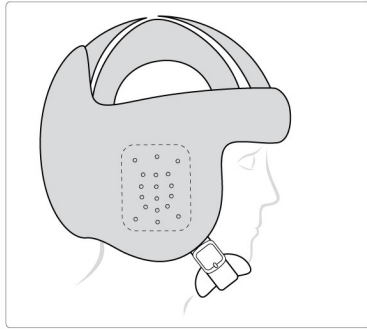


**...with alternative top
and ear protection**

Article number: 5791-...

**Chin protection for
Starlight® Protect-Evo**

Article number: 5555-01

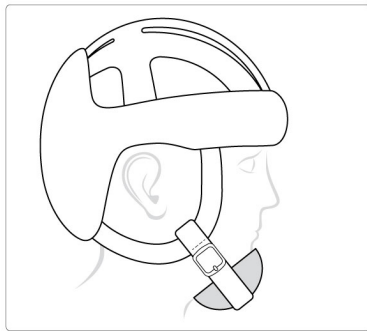

**Chin protection for
Starlight® Protect-Evo**

Article number: 5555-01



Strap system Biothane

Article number 5182-01BT

Particular robust
moisture-repellent and
easy to disinfect. Without
fastener. Type of fastener
freely selectable, except
fixlock fastener with
safety lock.




*Thorn buckle
(standard)*



*Strengthened
thorn
buckle*



*Velcro
fastener*

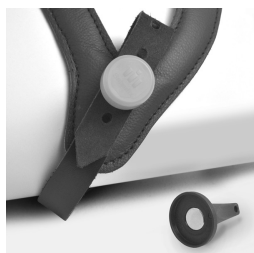

*Magnetic
fastener
including
key*



*Fixlock
fastener*



*Fixlock
fastener
with safety
lock*



Fasteners




*Thorn
buckle
(stand
ard)*



*Streng
thened
thorn
buckle*



*Velcro
fasten
er*

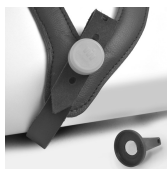

*Magne
tic
fasten
er
includi
ng key*



*Fixlock
fasten
er*




*Fixlock
fasten
er with
safety
lock*




*Fidlock
magne*


tic
snap
lock

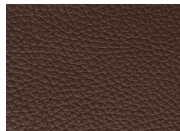


**Braun
(brown)**

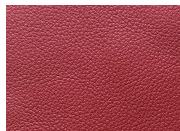

**1 - Rot
(red)**

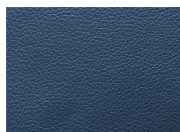

**2 -
Blau
(blue)**

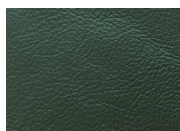


**4 -
Grün
(green
)**




**5 -
Schwa
rz
(black)**


Note! The colours shown are not binding!*
For more information please click here. → **Standard colours**

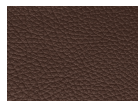


**Bra
un
(bro
wn)**

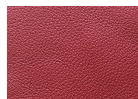


**1 -
Rot
(red
)**

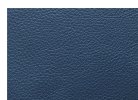

2 -
Bla
u
(blu
e)

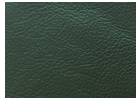

4 -
Grü
n
(gre
en)

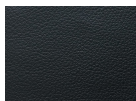

5 -
Sch
war
z
(bla
ck)

**Note! The colours shown are not
binding!***

For more information please click here. →

Please note:

When ordering please use the current ATO FORM

Starlight® - Dimension sheet!

Does your idea of a head protection look different?

You are welcome to send us your idea using the dimension sheet provided or communicate to us in a personal conversation.

*For technical reasons (including the calibration of the screen), it must be assumed that the original colours can only be represented approximately. Furthermore, deviations in colour nuances may occur due to different lighting conditions during product photography. Colour deviations are also possible during the finishing of the leather.

We therefore accept no liability for the colour of the product. Please also note our General Terms and Conditions (only in german).